
medivon.
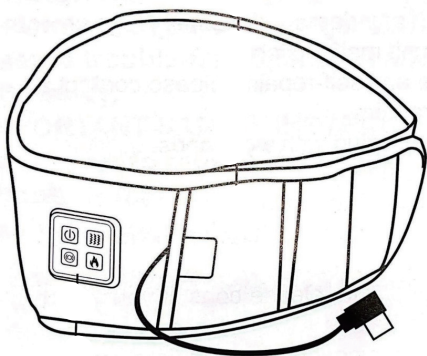
## User manual

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## Device manual

The device is ready for immediate use after purchase. Hold the power button for about 2 seconds to start the device.

After pressing the heating icon, the device will start warming up mode. Pressing the icon cyclically changes between the intensity (low - medium - high).

After pressing the massage icon, the device will start massaging the selected section. Pressing the icon cyclically changes between the intensity (low - medium - high).

The air massage icon starts a massage with an air cushion that gives deep relaxation and relaxation. Here are three cyclic modes ( low - med- high )

The device has a built-in 15-minute timer, after which it will turn off automatically. It is not recommended to massage one area of the body for more than 15 minutes to avoid injury and strain.

Note: Do not use the belt while lying down as this may cause personal injury or damage to the device

After the massage is completed, the device can be turned off by pressing the power button again.

Power source: the massager is powered by a classic USB-C cable (5V-2A). We recommend using a power bank to use the belt anywhere. You can also use a household power supply or a car socket. For proper operation of the device, a voltage of 2A is required, if it is lower, the massager will not allow you to run more than 1 function at the same time.

The strength of the massage should be adjusted to your own preferences. The device must not be used by persons with injuries or who have undergone surgery in the last 3 months. Do not use the device in high temperatures and damp rooms. Absolutely keep away from water. Also, people with hypertension, hypotension or severe strokes cannot use the device. People with skin problems, with a pacemaker or pregnant women should consult a doctor before use.

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Electrical and electronic devices as well as batteries and accumulators must not be mixed with other household waste. This is indicated by the symbol of the crossed out bin placed on the packaging. Worn-out electric and electronic equipment as well as batteries and accumulators should be delivered to local separate waste collection points or to the seller, and for details inquire in your commune. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage to them may result in damage to health or environmental pollution during further disposal or recycling of the equipment.

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The manufacturer hereby declares that the Medivon Orion device complies with Directive 2014/30 / EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility, and Directive 2011/65 / EU of the European Parliament and of the Council of 8 June 2011 on restriction of the use of certain hazardous substances in electrical and electronic equipment. The full text of the declaration of conformity is available at the internet address: [www.medivon.pl/deklaracje](http://www.medivon.pl/deklaracje)

Importer: Armare S.A. | ul Grzybowska 87, Warszawa, 00-844, Polska  
Model: Medivon Orion