
# BEAUTIFLY

SereneEye

PRO

series

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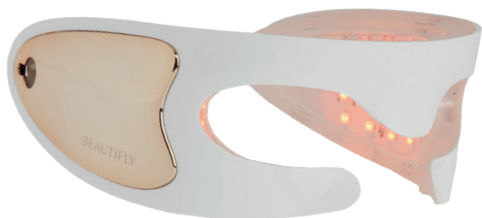- [www.beautifly.eu](http://www.beautifly.eu)

- [www.youtube.com/beautiflyeu](https://www.youtube.com/beautiflyeu)

Material:	ABS + silicone
Charging way:	USB-C charging
Charging time:	approximately 24 hours
Battery Capacity:	350 mAh
Tension:	3,7V
Output power:	1.2 W

#### Security:

Before using the product, read this user manual and follow the recommendations contained therein when using the device.

#### Device description [page 2-f1]:

On/off switch (1) Elastic band clasp (2) USB socket (3)

#### Light sensitivity test:

Before using the device for the first time, check whether you are allergic to light: place it on your forearm for 3 minutes. If skin becomes red or itchy, stop use and consult a dermatologist. Always clean your face before use to avoid reactions to cosmetics. If baking, stop using the appliance immediately.

#### How to use:

1. Cleanse your face and dry it.
2. Fasten and adjust the elastic headband.
3. Press the power button to turn on the device.
4. Press the mode button to select the desired mode (light therapy, micro vibration, light therapy + micro vibration)
6. After three minutes, the device will turn off automatically. Remove the mask and clean it
7. Apply a serum or moisturizing cream immediately after the treatment

#### Charging method:

The product should not be used while charging. (Note: The device will not work and buttons will not respond while charging.) It is recommended to fully charge the device before first use. Connect the USB-C cable to your device, then connect the USB plug to the adapter. When the power is turned on, charging will start. The device is fully charged when the indicator light turns from red to green. A fully charged device lasts approximately 30 minutes for 10 treatments.

#### Maintenance and care of the device:

Store this product in a dry place at room temperature and avoid exposure to dust, direct sunlight or humid conditions. Clean the mask before and after use.

#### Precautions:

Please keep the product away from small children and animals to avoid treating it as a toy. Handle the product carefully when using. This device is designed as a home electronic product and is not suitable for diagnosis or treatment. It should not be used by pregnant women, preparing for pregnancy or breastfeeding, people sensitive to light or taking photosensitizing drugs, with genetically determined skin diseases or conditions such as lupus, photosensitive eczema, or albinism. Do not use while charging. Do not use the device if it is damaged, Do not repair, disassemble or modify the device yourself. Any irregularities should be assessed by the service. Do not drop or hit the device to avoid malfunction. After each treatment, wipe the product with a soft damp cloth. Do not wash the device under water or submerge it. Do not use corrosive agents to clean the device. Store it in a cool, dry place.

#### More information about the company and products [page 2 - f2]

- [www.beautifly.eu](http://www.beautifly.eu)
- [www.youtube.com/beautiflyeu](https://www.youtube.com/beautiflyeu)

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**[EN] Rules for handling electrical and electronic equipment:** Waste electrical and electronic equipment as well as batteries and accumulators cannot be mixed with other waste. This is indicated by the crossed out bin symbol on the product packaging. Electrical and electronic equipment as well as batteries and accumulators may contain hazardous substances, mixtures and components harmful to the environment and human health. Improper handling or damage may result in damage to health or environmental pollution. Proper household handling of waste electrical and electronic equipment and batteries and accumulators positively contributes to the reuse and recovery, including recycling, of waste equipment. Waste electrical and electronic equipment as well as batteries and accumulators can be returned to local selective waste collection points. The distributor is obliged to collect used equipment from households free of charge at the point of sale, provided that the used equipment is of the same type and performed the same functions as the equipment sold. When delivering equipment intended for households to the buyer, the distributor is obliged to collect used equipment from households free of charge at the place of delivery of the equipment, provided that the used equipment is of the same type and performed the same functions as the delivered equipment. A distributor running a retail unit with a sales area within the meaning of Art. 2 point 19 of the Act of March 27, 2003 on spatial planning and development of at least 400 m<sup>2</sup> devoted to the sale of equipment intended for households, is obliged to accept waste equipment from households free of charge in this unit or in its immediate vicinity, none of the external dimensions exceeds 25 cm, without the need to purchase new household equipment. The product was introduced to the market on August 13, 2005

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**[EN] Summary declaration of conformity:** The manufacturer declares, under its sole responsibility, that the SereneEye PRO product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonization of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Directive Commission delegated act (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of substances subject to restriction and Directive (EU) 2017/2102 of the European Parliament and of the Council of November 15, 2017 amending Directive 2011/65/EU on the restriction of the use of certain hazardous substances in electrical and electronic equipment. The full content of the declaration of conformity is available at [www.beautifly.eu/deklaracje](http://www.beautifly.eu/deklaracje)

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