# BEAUTIFLY

SkinRenew

P R O


series

DermaPen  
DermaPen

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f1

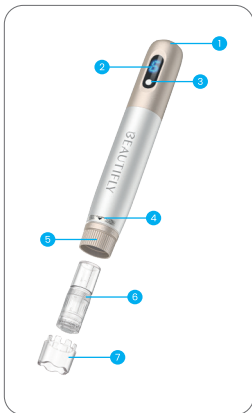f2

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Input:	100 - 240V 50/60Hz
Output:	5V $\bar{\bar{}}$ 2A
Number of modes:	6
Included:	Dermapen device, needle cartridge, charger.

#### Security:

Before using the product, read this user manual and follow the recommendations contained therein when using the device.

#### Device description [page 2 - f1]:

1. Charging input.
2. Intensity display.
3. On/off/speed button.
4. Scale point.
5. Scale adjustment ring.
6. Part of the cartridge with needles.
7. Protective cover.

#### Landing:

1. This device is powered by a built-in battery. Can work wirelessly. Charging the device takes about 2 hours.
2. The power indicator remains red while charging.
3. The power indicator will turn green when charging is completed.
4. The device speed number flashes slowly when the battery level is 10%. Once fully discharged, the device will turn off automatically and the speed level will flash on the display for one minute.
5. The device shows "L" symbol when restarting after automatic shutdown.
6. The device shows "0" symbol while charging.
7. The device shows "F" symbol when charging is complete.

#### Preparation:

1. Remove the needle cartridge from the sterilized packaging.
2. Align the round socket with the machine, press the cartridge horizontally until it stops.
3. Press the power button for 2-3 seconds to turn on or turn off the device.

#### Needle length adjustment:

1. Adjust the needle depth by turning the scale scale adjustment ring.
2. Turn it counterclockwise, the needle depth will be shorter.
3. Turn it clockwise, the needle depth will be longer
4. Each time you press the speed button, its mode changes.

#### Recommendations:

1. Used needle cartridges should be placed in the garbage container.
2. If the needle cartridge is accidentally contaminated during installation (dropped on the floor, package open/damaged), discard it and purchase a new needle cartridge.

#### Contraindications:

- active inflammation and purulent processes during the procedure
- pregnancy and breastfeeding
- unregulated diabetes

- unregulated thyroid diseases
- tendency to develop keloid scars
- allergy to any of the ingredients of the substances used for the treatment
- in people with vascular skin
- use of drugs that reduce blood clotting
- active cancer diseases
- use of anticoagulant drugs (e.g. aspirin)
- fever
- if Botox has been administered to the treated area in the last three months

Warnings and recommendations:

1. Please check the device and charger before use.
2. Make sure the device is turned off when not in use.
3. Keep the device away from direct sunlight and high temperature or moisture.
4. Keep the device out of the reach of children.
5. Turn off the device after use or before cleaning and maintenance.
6. Do not use the device on the eyelids, thyroid area, lips, ears, breasts or genital area.
7. Do not repair, disassemble or modify the controller and body. This may cause a malfunction.
8. When using a device that has not been used for a long time, check whether it is working properly and safely.

More information about the company and products [page 2 - f2]:

- [www.beautify.eu](http://www.beautify.eu)
- [www.youtube.com/beautifyeu](https://www.youtube.com/beautifyeu)

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**[EN] Rules for handling electrical and electronic equipment:** Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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**[EN] Summary declaration of conformity:** The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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