# BEAUTIFLY

MicroGlow

PRO


series

DermaPen  
DermaPen

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

f1

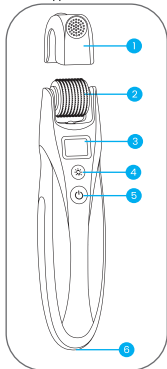f2

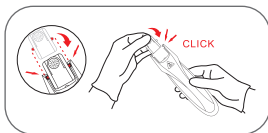f3

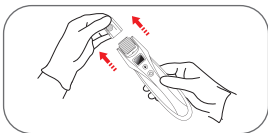f4

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Input:	100 - 240V 50/60 Hz
Output:	5V $\approx$ 2A
Number of modes:	5
Included:	Dermapen device, two massage rollers, charger.

#### Security:

Before using the product, read this user manual and follow the recommendations contained therein when using the device.

#### Device description [page 2 - f1]:

1. Protective cover.
2. Microcurrent roller.
3. Display screen.
4. Backlight button/vibrate button.
5. ON/Off button, change the intensity.
6. Power port.

The device combines microneedling technology with microcurrents.

**Advantages of microneedling:** Improves the appearance of the skin, reduces the signs of aging. Stimulates the skin using microneedles to produce collagen and elastin. The micro-channels created by the device on the skin's surface also aid in the absorption of serum.

**Advantages of microcurrent:** Increases cell activity, stimulates muscle movement and restores elasticity, accelerates blood circulation in capillaries.

#### How to use [page 2 - f2, f3]:

1. Remove the needle cartridge from its packaging.
2. Hold the device with the word "SKIN" facing up. Then install the shaft with the protective cap on the handle connector. (Before installing the roller on the handle, it must be sterilized.)
3. Remove the protective cap from the roller. The device is now ready for use.
4. Press and hold the power button for 3 seconds to turn on or turn off your device.
5. Use the appropriate intensity level depending on the treatment and skin type (there are 5 different levels). Change the intensity of the microcurrent by pressing the intensity button.
6. Press the light button to turn the light and vibration on or off.
  - blue light - weak vibrations, works well in cleansing the skin and helps in the fight against acne.
  - red light - reduces inflammation and stimulates collagen production.

#### Landing:

1. This device is powered by a built-in battery. Can work wirelessly. Charging the device takes about 2 hours.
2. The device speed number flashes slowly when the battery level is 10%. Once fully discharged, the device will turn off automatically and the speed level will flash on the display for one minute.
3. The device shows "L" symbol when restarting after automatic shutdown.
4. The device shows "O" symbol while charging.
5. The device shows "F" symbol when charging is complete.

**Contraindications:**

- active inflammation and purulent processes during the procedure
- pregnancy and breastfeeding
- unregulated diabetes
- unregulated thyroid diseases
- tendency to develop keloid scars
- allergy to any of the ingredients of the substances used for the treatment
- in people with vascular skin
- use of drugs that reduce blood clotting
- active cancer diseases
- use of anticoagulant drugs (e.g. aspirin)
- fever
- if Botox has been administered to the treated area in the last three months

**Warnings and recommendations:**

1. Please check the device and charger before use.
2. Make sure the device is turned off when not in use.
3. Keep the device away from direct sunlight and high temperature or moisture.
4. Keep the device out of the reach of children.
5. Turn off the device after use or before cleaning and maintenance.
6. Do not use the device on the eyelids, thyroid area, lips, ears, breasts or genital area.
7. Do not repair, disassemble or modify the controller and body. This may cause a malfunction.
8. When using a device that has not been used for a long time, check whether it is working properly and safely.

More information about the company and products [page 2 - f4]:

- [www.beautifly.eu](http://www.beautifly.eu)
- [www.youtube.com/beautiflyeu](https://www.youtube.com/beautiflyeu)

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**[EN] Rules for handling electrical and electronic equipment:** Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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**[EN] Summary declaration of conformity:** The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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