BEAUTIFLY

Aura

PRO

series

Maska LED
LED Mask

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beautifly.eu

@beautifly.eu

beautifly.eu


/beautiflyeu

@beautifly.eu

f1

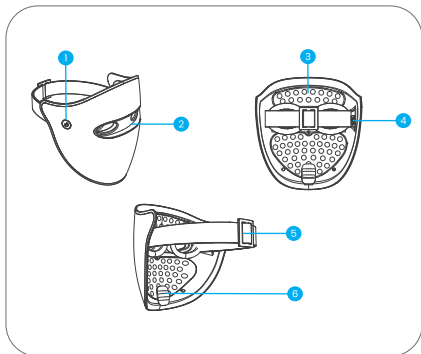f2

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- www.beautifly.eu
- www.youtube.com/beautiflyeu

Rated voltage:	DC3.7V
Working current:	Soothing mode<600mA, Brightening mode<900mA, Anti-aging mode 1300mA
Battery Capacity:	1000 mAh
Charging time:	<3H
Smart Timer:	10 minutes
Standby time:	Soothing mode>80 minutes, Brightening mode>60 minutes, Anti-aging mode 40 minutes
Photon wavelength:	Deep red light: 660±5nm, red light: 630±5nm, infrared light: 850±5nm, blue light: 470±5nm
Irradiation:	Soothing mode 55±5mW/cm ² , Brightening mode 55±5mW/cm ² , Anti-aging mode 75±5mW/cm ²

Security:

Before using the product, familiarize yourself with the operating instructions and follow the instructions when using the device.

Device description [page 2-f1]:

1. ON/OFF button, lamp control.
2. The window will fall off your eyes.
3. Photon diode.
4. Charging port.
5. Adjustable headband.
6. Chin frame.

Functions:

1. Soothing mode: press button once, start indicator light flashes white, beep sounds, voice prompt sent "Soothing mode". The blue light will remain on and the device will turn off automatically after 10 minutes.
2. Brightening mode: apply the power button, the start indicator light flashes white, it will beep, output voice prompt "Brightening mode". Red light and light is on and the device turns off automatically after 10 devices.
3. Anti-aging mode: press the button three times, the indicator light will flash white, the beep will sound, and the voice prompt will display "Anti-aging mode". Red light, deep red color and infrared light, on, the device turns on automatically after 10 minutes.

Note: If no operation is performed, the device will automatically turn off after 10 operations. Duration does not accumulate when switching modes while running.

How to use:

1. Cleanse your face and dry it.
2. Put on the mask and freely choose the skin protection mode.
 - soothing mode is intended for skin with unstable water-fat balance, such as or moisture mobility
 - the brightening mode is intended for dull and dull complexions
 - anti-aging mode is on device with fine lines, wrinkles or available
3. After 10 activities, the mask turns off automatically.

During use, adjust the chin frame to the shape of your face and adjust the bands to fit.

Battery:

1. Turn off/press to turn off. It will beep twice, the light will turn off and the device will stop working.
2. When the battery power level drops below 30%, the power button flashes red. The device can still be used for 10 minutes. When the battery level is critically low, the device cannot be turned on. Application/indicator light flashes red, produces two short beeps.
3. When disconnected, the disconnect/off button/light remains white and turns off when disconnected.

Cleaning:

After use, clean the safety device and basic cloth. Do not spray the mask in water, feed it dry.

Measures to be taken:

1. Not applicable if in the next case the effects or effects that will occur.
2. Do not use the device if you have inflamed, scarred, cracked skin or are allergic to electromagnetic fields.
3. It is not intended for use by pregnant women.
4. Do not look directly at the beauty light before confirming that the product is ready to use.
5. Do not create the device while it is running.
6. Do not use the device while charging.
7. Do not disassemble the device.
8. Do not consider the device if the silicone in the eye shield is a light source.
9. Connection that the cable and power supply are in good condition to avoid the risk of electric shock.
10. If the product is hit or dropped, please check for any wear before use.
11. Stay away from high temperature environment.
12. stay away from the presenter.
13. Stay out of reach of children.
14. If the device has not been used for a further period of time, it is recommended to fully charge it before use.

More information about the company and producers [page 2 - f2]:

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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Importer: Beautify Sp. z o.o.

Plac Bankowy 2, 00-095 Warszawa, Polska

Model: Aura