BEAUTIFLY

Hydra Peel

PRO

series


Hydrabrasion device
User manual

 /beautifly.eu  @beautifly.eu  beautifly.eu

 /beautiflyeu  @beautifly.eu

f1

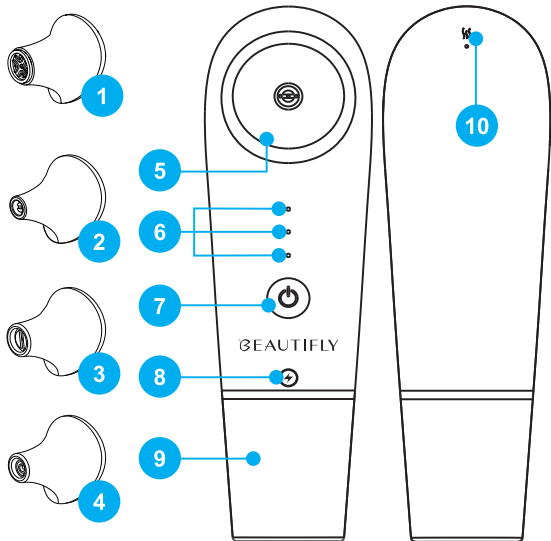f2

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- www.beautifly.eu

- www.youtube.com/beautiflyeu

Product color:	white, rose gold
Device weight:	approx. 166 g
Device dimensions:	188 x 55 x 75 mm
Rated voltage:	3.7V
Supply voltage:	5V — 1A
Rated power:	1W
Intensity levels:	3 levels
Temp. device storage:	0oC – 45oC
Working environment temperature:	5oC – 35oC
Working environment humidity:	25% – 75%

Security:

Before using the product, read this user manual and follow the recommendations contained therein when using the device.

Device description [page 2 - f1]:

1. Silicone cleaning tip.
2. Small cleaning tip.
3. Large cleaning tip.
4. Spiral cleaning tip.
5. Cleaning tip.
6. Intensity level indicator.
7. Power button (On/Off switch).
8. Charging socket.
9. Container.
10. Ventilation.

Hydrogen purification (hydrabrasion)

Hydrabrasion is a gentle treatment that cleanses, moisturizes and exfoliates the skin. It is definitely more delicate than microdermabrasion, so there are no contraindications to its use. Even people with sensitive, mature, vascular or acne-prone skin can benefit from this type of treatments. In order to achieve maximum cleansing effects, use special tips for treatments on specific areas of the face. Thanks to them, you will not only cleanse the skin even in hard-to-reach places, but also easily achieve the expected results.

Device operation:

1. Press and hold the power button for 3 seconds to turn your device on or off.
2. Briefly press the power button to increase or decrease the intensity level.
3. The first intensity level turns on one green light, the second intensity level turns on two green lights, and the third intensity level turns on three green lights.
4. While charging, the charging indicator lights red. When the battery is full, the red light turns off and the green light turns on.
5. It takes approximately 3.5 hours to fully charge the device.

Precautions during use:

1. Do not place the device in an environment with a strong magnetic field.
2. Be careful not to tilt the device more than 45°C when using, otherwise water will spray out of the exhaust port.
3. This device can only use clean water. Do not use other liquids that are not suitable for this device.
4. Do not use the device while it is charging.
5. The inside of the water tank can be disinfected and sterilized.

6. This device should not be soaked or immersed in water for cleaning. Be careful when using the device and ensure that no water penetrates the device when using or cleaning it.
7. After each use of the device, empty the water tank of used water and then dry the device thoroughly.

How to use the device:

Remove makeup from your face using a makeup remover. Install the filter cap on the bottom of the device and add 2/3 clean water to the water tank. Attach the appropriate cleaning tip. There are different sizes of cleaning tips. Depending on your needs, you can choose the right cleaning tip to use in different positions. After cleansing your face, wash it with tonic and apply a moisturizer to your skin.

Description of the functions of the cleaning tips:

Spiral cleaning tip (small) - Suction center (commonly used), Function: dead skin, cuticle, pore cleansing (universal type), Location: small area Area: T zone, nose part.

Cleaning tip (large) - Large suction (commonly used), Function: removal of blackheads, dead skin, calluses, cleansing pores (universal type),

Location: large surface of the face, Area: cheeks, forehead.

Cleaning tip (small) - Enhanced suction power (persistent contamination),

Function: blackheads, pimples, cuticles, pore cleansing (pre-blackheads) (universal type), Location: small area, Area: T zone, nose part.

Silicone cleaning tip (large) - Enhanced suction power (persistent dirt), Function: acne, calluses, cleansing pores (preferred blackheads) (universal type), Location: large surface of the face, Area: cheeks, forehead.

More information about the company and products [page 2 - f2]:

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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Shortened warranty terms

This Guarantee is granted by a Polish company, Beautifly Sp. z o.o. based in Warsaw. Plac Bankowy 2, 00-095 Warsaw, Poland. You can find us in the official Polish register of enterprises under the KRS number: 0000887306, NIP: 5252854218 and REGON: 388382021 We guarantee that the equipment you have purchased from us (hereinafter referred to as the "Product") will work properly. If something goes wrong for a reason inherent in the Product, this warranty protects you. The warranty starts from the day you receive the Product and lasts:

24 months, if you are a consumer or a natural person running a business and the purchase of the Product is not of a professional nature for you. Otherwise, the warranty period is 12 months. If the law in your area requires a longer warranty, the warranty period will be extended to meet local law requirements.

In addition, if the Product does not comply with the contract, you have the right to legal protection from and at the expense of the seller. This warranty does not affect these rights - it does not limit, suspend or exclude them.

Detailed warranty conditions, how to submit complaints, what happens when you find a defect, what are exceptions to the warranty, when the warranty expires or is extended, and what are the rules regarding the processing of personal data, you can find on our website: www.beautifly.eu/warranty-beautifly