S·E·Y·S·S·O

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BEFORE YOU BEGIN


Thank you for choosing the SEYSSO DINOSAUR TRAVEL water flosser. Please carefully read this manual before using it for the first time to ensure optimal performance and safety. Remember, this product is designed for children over the age of four and should only be used under adult supervision. Keep this manual for future reference.

PACKAGE CONTENTS

- 1 SEYSSO Dinosaur Travel device
- 1 USB cable
- 2 SEYSSO Dinosaur Standard tips
- 1 user manual

DEVICE DESCRIPTION

1. Power button
2. Mode selection button
3. LED indicator
4. Water tank
5. Charging port

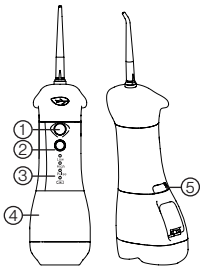USAGE SAFETY

1. Use the device as intended. Do not attempt to modify or disassemble the device.
2. Store the device away from moisture and high temperatures.
3. Do not use the device if the cord or plug is damaged.
4. If smoke, excessive heat, or strange odors occur, immediately turn off the device and contact customer service.

WARNINGS! READ BEFORE USE! KEEP THIS MANUAL

1. For household use only. 
2. The flosser is equipped with a waterproof handle/body.
3. It is safe and designed for bathroom use.
4. The device has an IPX7 protection rating, meaning it is resistant to accidental submersion in water up to 1 meter deep for up to 30 minutes.
5. Reminder: the flosser is not designed for underwater operation!
6. The product is not intended for self-repair.
7. Unwind the cable before use.
8. In case of failure, do not disassemble or replace any components, as there is a risk of electric shock. Contact the manufacturer's service department.
9. A damaged power cord must be replaced with a new one. Contact the manufacturer's service department.
10. If the flosser is dropped, replace the tip before the next use, even if no damage is visible.

MAINTENANCE AND CLEANING

1. Remove the tip from the flosser and rinse it under running water after cleaning the interdental spaces.
2. Do not use water hotter than 40°C or acidic cleaning solutions.
3. Wipe the flosser dry before connecting it to a charger. Ensure the flosser is unplugged before cleaning. Use a dry cloth to wipe the charger surface.
4. Follow basic precautions when using any electrical device, especially around children.

CHARGING THE DEVICE

1. Connect the USB cable to the charging port on the device.
2. Plug the other end of the cable into a USB power adapter.
3. The LED indicator will light up to indicate charging has started.
4. Full charging takes approximately 4-6 hours.
5. The flosser is charged using a USB C-USB A 5V cable.
6. Each subsequent charging session should last around 4-6 hours.
7. Disconnect the flosser from the cable while cleaning teeth.
8. There is no risk of overcharging the battery.
9. To extend battery life, charge the flosser when the low battery indicator (red light) illuminates.
10. Once fully charged, the battery indicator will turn solid green.

DANGER WARNING

To reduce the risk of electric shock:

1. Always disconnect the device from the power source after use.
2. Do not use during bathing.
3. For more detailed instructions on storing and using the device, as well as the included accessories, refer to the manufacturer.

To reduce the risk of fire, short circuit, or injury:

1. Disconnect the charger from the power source after charging.
2. This product is not a toy. Use extreme caution if the flosser is used by or near a child.
3. Use only accessories recommended by the manufacturer.
4. Do not use a damaged charging cable. If the flosser does not operate correctly after dropping, exposure to the elements, or submersion in water, return it to the point of purchase.
5. Do not stretch, pinch, or wrap the charging cord around sharp edges. Do not use the flosser while charging. Keep device components away from heated surfaces.
6. Do not plug the cord in with wet hands.
7. Never use the product if you are drowsy or unable to concentrate.
8. Do not use the device outdoors or in environments exposed to aerosols or pure oxygen.

ENVIRONMENTAL PROTECTION

Do not dispose of the SEYSSO Dinosaur water flosser with household waste. Take it to an electronic recycling and disposal point.

DEVICE USE

1. Press the power button to turn on the device.
2. Toggle modes by repeatedly pressing the mode selection button until the desired mode is selected.
3. Turn off the device after use by pressing the power button again.
4. Regularly clean the device to ensure long-term performance.

TROUBLESHOOTING

- **The device does not turn on:** Ensure it is fully charged.
- **Charging takes longer than usual:** Check if the charging cable and power adapter are functioning correctly.
- **The LED indicator does not light up during charging:** This may indicate an issue with the charging port or cable.

STORAGE

1. Store the device in a dry place, away from direct sunlight.
2. Avoid extreme temperatures, which can affect device performance.
3. If you do not plan to use the device for an extended period, charge it fully before storing.

CONTACT

For questions or technical issues, please contact our customer service department:

Customer service: seysso@seysso.com

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Website: www.seysso.pl

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Model	SEYSSO DINOSAUR TRAVEL
Input Voltage	5V  1A
Charging Time	Approximately 4-6 hours
Battery Life	Up to 90 minutes of use after a full charge
Dimensions	100 mm x 50 mm x 20 mm
Weight	150 g
Charging Method	USB-DC 5V cable
Rated Input Power	3W
Battery	Lithium-ion
Battery Capacity	1400 mAh
Low Battery Indicator	LED light
Waterproof Rating	IPX7

Important: Use of an Appropriate Power Adapter

Ensure that the power adapter used has output parameters consistent with those specified on your flosser's label. Using a power adapter with incompatible parameters may result in improper device operation and is strictly prohibited. Failure to comply with this guideline may cause device damage and void the warranty. Take care of your device by using only the recommended accessories and power adapters.

WARRANTY

SEYSSO provides a two-year warranty from the date of purchase. The warranty covers damages related to manufacturing defects or improper product craftsmanship. Retain the proof of purchase for the duration of the warranty. During the warranty period, you may file a claim if the product has been used correctly and as intended. THE WARRANTY DOES NOT COVER:

1. Surface damages resulting from normal use.
2. Damages resulting from:
 - a) Accidental dropping or impact, improper use, abuse, commercial use, negligence, or modification of the device.
 - b) Improper maintenance or lack of care for the product.
 - c) Connection to an incorrect power source.
 - d) Repairs or attempted repairs conducted by a party unauthorized by SEYSSO.

Any attempt to open or disassemble the device, its parts, or accessories will void the warranty. Warranty claims must be supported by evidence and must be submitted within the warranty period.

The material covering the handle/body is delicate; therefore, care should be taken during use. Scratches or abrasions on the surface of the handle will be considered surface damage resulting from normal use and will not be covered by the warranty.

The warranty applies to components that affect the functionality of the flosser.

This manual applies to models with the following numbers: SK.WF.010.COR.1

MANUFACTURER:
WM Brands sp. z o. o.
ul. Dziupli 12
02-454 Warsaw, Poland
Made in P.R.C.

CONTACT:
Customer service: seyso@seyso.com
Warranty claims: reklamacje@seyso.com
Website: www.seyso.pl

Water Resistance (Flosser Handle): IPX7. Protection against submersion in water up to a depth of 1 meter for 30 minutes.

The flossers are equipped with a built-in lithium-ion battery.

The product complies with the requirements of the European Union (EU) „New Approach” directives.

The packaging is recyclable.

The manual is recyclable.

This device is marked in accordance with the Polish Act on Waste Electrical and Electronic Equipment and Directive WEEE 2012/19/EU with the symbol of a crossed-out wheeled bin. This marking indicates that the device, after its period of use, must not be disposed of together with other household waste.

The user is obliged to return it to entities collecting waste electrical and electronic equipment. Collection entities, including local collection points, stores, and municipal units, have established a system that allows for the proper disposal of this equipment. Proper handling of waste electrical and electronic equipment helps avoid harmful consequences for human health and the environment caused by the presence of hazardous components and improper storage and processing of such equipment.

This device has been designed to limit the presence of certain substances that may negatively impact the environment during its use and after its disposal.

For household use only