S · E · Y · S · S · O

MODEL: SH.ST.030.NVB.1

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- 1 brush head
- 1 brush handle
- 1 USB cable

Connecting and charging

- The brush is charged via a 5V USB-C cable.
- The first full charge takes about 6 hours and lasts up to 10 months (in Sensitive mode).
- Each subsequent charge should last about 3.5-6 hours.
- During charging, the brush must be disconnected from the cable if you want to brush your teeth. It is not possible to brush your teeth while the brush is connected to the charger. After brushing your teeth, you can reconnect the brush to continue charging.
- There is no risk of overcharging the battery.
- To prolong battery life, it is recommended to charge the brush when it signals a low battery level ($\leq 20\%$).
- Once fully charged, the display will show „100%“.

Using the sonic toothbrush

Place the head on the handle. Then apply the appropriate amount of toothpaste to the brush bristles. Before turning the brush on, insert the head into your mouth. Turn the brush on. After the appropriate time, depending on the selected program, the device will automatically turn off and be ready for reuse. The toothbrush, during tooth brushing, informs you to change the brushing zone every 30 seconds on the 2-minute mode. The brush moves at different speeds, depending on the selected mode:

Sensitive: 56,000 sonic movements per minute

Clean: 74,000 sonic movements per minute

White: 82,000 sonic movements per minute

DIY: 96,000 sonic movements per minute

In each mode, you can choose the intensity level from 1 to 5.

NOTE: The brush bristles will start working only when the brush is properly pressed against your teeth.


Usage notes

- After use, rinse the brush head with clean water and let it dry naturally.
- Replace the brush head only with the original SEYSSO Hypertech Infinity Brush Head (replaceable head). The brush functions in all modes only with this original head. Using other heads may cause the brush to work improperly and it will function only in DIY mode.
- You can rinse the brush directly with water, but do not use water above 40°C or acidic cleaning solutions.
- Fully discharge the brush at least once a year.
- If the charging cable is damaged, please contact customer service.
- If the product falls, the handle or brush head may be damaged. Do not use it after damage.

- The body or charging port must be kept away from heat sources and should not be used in high temperatures.
- Do not disassemble or attempt to repair the device yourself. Contact post-sales service if repairs are needed.
- If you do not plan to use the device for an extended period, rinse it with water and store it in a cool, dry place, away from direct sunlight.

Product part names

- A. Brush head
- B. Battery charge indicator
- C. USB charging port
- D. Handle
- E. LCD display
- F. Function selection button
- G. Power and function confirmation button
- H. Pressure indicator

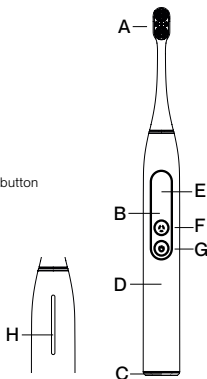Usage instructions

Step 1: Click to select the brushing mode.


Step 2: Enter the brushing mode, and the handle display will show the selected program and its brushing strength.

Step 3: After brushing, the handle display will show the battery and brushing status.

Brushing guide

Step 1: Place the brush head onto the peg, do not rotate it.

Note: It is normal for there to be a gap between the brush head and the handle.

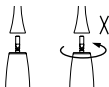Step 2: Lightly wet the brush bristles and apply toothpaste.

Step 3: Place the brush head at a 45° angle to the tooth arches and gently press so the bristles touch the gum line or just below the gum line.


Note: The center of the brush head should always be in contact with the teeth.

Step 4: Click the power button to turn on the brush.

Step 5: To clean the inner surfaces of the front teeth, hold the brush vertically.

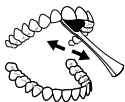Step 6: This brush divides the oral cavity into 16 zones, including internal, external, and chewing surfaces. It can recognize the zones being brushed and remind you of missed areas to ensure comprehensive cleaning of the oral cavity.

1. Internal incisors
2. External incisors
3. External molars
4. Internal molars
5. Chewing surfaces of molars

Technical features

a. Mode selection

- Sensitive: Suitable for sensitive teeth.
- Clean: Daily oral hygiene.
- White: Effective stain removal from teeth.
- DIY: Your favorite mode that cleans the oral cavity.

For each mode, you can set the working time: 1 minute, 2 minutes, 3 minutes, unlimited.

b. Intensity selection

You can select different intensity levels according to your preferences.

c. Pressure feedback

The product can detect the brushing force using a real-time pressure monitoring system and cooperate with a gum bleeding detection algorithm. It provides feedback through a light indicator on the back of the handle, handle vibration, and an audible warning to intelligently remind you to control the pressure.

- Blue: Pressure too weak.
- Green: Pressure appropriate.
- Red: Pressure too strong.

d. Brush head replacement reminder

- The brush head replacement reminder may appear on the display.

The toothbrush will remind you to replace the brush head after 180 usage cycles. Once this reminder is displayed, the usage counter will be reset to zero, and the toothbrush will start recording each use again.

Usage recording details:

- Each time the toothbrush is used for more than 30 seconds, it is counted as one usage.
- Uses shorter than 29 seconds will not be recorded as a usage.

e. Improper brush head warning

- During the installation of the brush head, if the device detects a non-original head, a corresponding warning will appear.

f. Shock absorption

- The product has soft rubber shock absorption between the motor position and the body, ensuring that the brush does not cause discomfort during use.

Battery status and charging

a. Battery status

- After turning it on, if the battery level is too low, the display will show:

- Power indicator below 20%, low battery reminder, which can be removed by pressing the mode button.
- If the battery is below 10%, the screen will automatically turn off after 5 seconds.

b. Charging status


- During the charging process, the display will show the charging status in real time.

c. Charging instructions

Step 1: Ensure the brush is dry, and clean it if wet.


Step 2: Connect the USB cable to the adapter. Adapter specifications: 5V 1-2 A. Then connect the adapter to the power outlet.

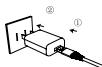Step 3: Insert the USB cable into the brush charging port; if charging is successful, the handle display will light up.

Step 4: Wait until the brush is fully charged.

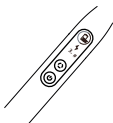WARNINGS! READ BEFORE USE!

KEEP THIS MANUAL

- For home use only.
- The toothbrush has a waterproof handle. It is safe and designed for bathroom use. The device has an IPX7 protection rating. This means the toothbrush is resistant to accidental immersion in water up to 1 meter deep for up to 30 minutes.
- Remember: the toothbrush is not intended for underwater operation!
- The product is not intended for self-repair.
- Straighten the cable before use.
- In case of malfunction, do not disassemble the device or replace any parts, as there is a risk of electric shock. Contact the manufacturer's service department.
- A damaged power cord must be replaced with a new one. Contact the manufacturer's service department.
- If the toothbrush is dropped, replace the head before using it again, even if damage is not visible.

Maintenance and cleaning

After brushing your teeth, rinse the toothbrush under running water. Pay special attention to the sealing areas: near the peg and around the button, where toothpaste can accumulate and promote the growth of microorganisms. Wipe the handle with a soft cloth.

- Remove the head from the handle and rinse it under running water.
- Do not use water hotter than 40°C or acidic cleaning solutions.
- Wipe the handle before connecting it to charging.
- Before cleaning the toothbrush, ensure it is unplugged from the power supply. Use a dry cloth to wipe the surface of the charger! When using any electrical appliance, especially in the presence of children, take basic safety precautions.

Attaching and removing the brush head

Before the first use of the toothbrush, apply a thin layer of lubricant to the metal peg of the handle. This will make it easier to attach and remove the heads. It is recommended to replace the heads every 3 months.

Danger warnings

To reduce the risk of electric shock:

- Always disconnect the device from the power source after use.
- Do not use while bathing.
- For more detailed instructions on storage and use of the device, as well as the included accessories, contact the manufacturer.

To reduce the risk of fire, short circuit, or injury:

- After charging is complete, unplug the charger from the power source.
- The product is not a toy. Use extreme caution if the toothbrush is used by a child or in the presence of a child.

- Use it only as intended. Use only accessories recommended by the manufacturer.
- Do not use the charging cable if it is damaged. If the toothbrush has fallen, been left outdoors, or submerged in water and does not work properly, return it to the place of purchase.
- Do not pull the charging cable, cut it with doors, or wrap it around sharp edges. Do not turn on the toothbrush while charging. Keep the device's components away from heated surfaces.
- Do not hold the plug with wet hands.
- Never use the product if you have trouble concentrating or feel sleepy.
- Do not use the device outdoors or in areas exposed to aerosols or pure oxygen.
- Before plugging the device into an outlet, familiarize yourself with its specifications.

It is not recommended for children under 3 years old to use the device. Children and people with reduced sensory, motor, or mental capabilities should use the device under the supervision of a guardian or be properly trained in its correct and safe use. Children should not play with the device. Children should not perform cleaning and maintenance of the device unsupervised. Do not use excessive force or quick movements while brushing your teeth. In the first few days of using the toothbrush, the gums may slightly bleed. Bleeding should stop on its own after a few days. If it persists or worsens, we recommend contacting a dentist. For sensitive gums and teeth, it is recommended to use the Soft mode.

Environmental protection

Do not dispose of the SEYSSO Infinity toothbrush with household waste. Take it to an electrical and electronic equipment recycling center.

Two-year warranty

SEYSSO provides a two-year warranty from the date of purchase. The warranty covers damage related to manufacturing defects or improper execution of the product. Keep the proof of purchase for the duration of the warranty. During the warranty period, you may file a claim if the product was used correctly and in accordance with its intended purpose.

THE WARRANTY DOES NOT COVER:

1. Surface damage resulting from normal use.
2. Damage resulting from:
 - a) Accidental dropping or impact, improper use, abuse or commercial use, negligence, or device modification,
 - b) Improper maintenance, negligence toward the product,
 - c) Incorrect connection to the power source,
 - d) Repairs or repair attempts by parties not authorized by SEYSSO. Any attempt to open or disassemble the device, its parts, or accessories will void the warranty. Warranty claims must be supported by evidence. The claim must be made within the warranty period. The material covering the handle is delicate, so it should be handled with care. Scratches or wear on the surface of the handle will be considered surface damage resulting from normal use.

Therefore, they will not be covered by the warranty.

The warranty covers components that affect the operation of the toothbrush.

Technical specifications

PRODUCT NAME	SEYSSO HYPERTECH INFINITY SONIC TOOTHBRUSH
Switch type	2 buttons
LCD display	LCD
Operating time	Up to 1200 minutes
Charging method	USB-C 5V cable charging
Charging time	Approximately 3.5–6 hours
Operating voltage	3.7V
Rated input power	≤5W
Battery	Lithium-ion battery
Battery specification	Lithium-ion /18650 /3000mAh

Important: Use the appropriate charger. It is recommended to use a charger with a standard charging current of 5V - 1A or 2A. Ensure that the charger you are using has output parameters matching the ones stated on your toothbrush label. Using a charger with incompatible parameters may result in improper device operation and is strictly prohibited. Failure to comply with this rule may result in device damage and loss of the warranty. Take care of your device by using only recommended accessories and chargers.

The instructions apply to model number: SH.ST.030.NVB.1

MANUFACTURER:

WM Brands sp. z o.o.

ul. Dziupli 12

02-454 Warsaw, Poland

Made in P.R.C.

CONTACT:

customer service: seysso@seysso.com

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www.seysso.pl

Water resistance (brush handle): IPX7. Protection against immersion in water up to a depth of 1 meter for 30 minutes.

The toothbrushes have a built-in lithium-ion battery.

The packaging is recyclable.

The product complies with the requirements of the European Union „New Approach” directives.

The instruction manual is recyclable.

This device is marked with the symbol of a crossed-out waste bin in accordance with the Polish Act on Waste Electrical and Electronic Equipment and the WEEE Directive 2012/19/EU.

This marking informs that the device, after the end of its service life, must not be disposed of with other household waste. The user is obligated to return it to an electrical and electronic waste collection facility. Collection facilities, including local collection points, stores, and municipal units, create a system enabling the return of this device. Proper handling of electrical and electronic waste helps prevent harmful effects on human health and the environment caused by the presence of hazardous components and improper storage and processing of such equipment.

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